



Green Button Go™



Lab Automation with Compliance Control

Digital records have been crucial to expediting information flow in the lab. If you work in a highly regulated lab, tracking and validating electronic records and signatures is required to comply with the Code of Federal Regulations Title 21 for the Food and Drug Administration, also known as 21 CFR Part 11.

Green Button Go™ Scheduling Software is available in a version that is certified 21 CFR Part 11 compliant, so you can maintain security and traceability in the lab while enjoying seamless laboratory automation. The compliant software ensures locked-in access, audit trails, electronic signature logs and other required documentation for using electronic records.

Control Access

- Set different permission levels for each user account
- Password protection options include password rules and expiration intervals

Validate Programs

- Only administrators can establish a validated program or make edits to it
- Lock down access to the program so users can't make changes
- System users can only run a validated program

Log Electronic Signatures

- Users can sign a log file or Sign at the end with an electronic signature
- Signatures on the log files are recorded

Record Audit Trails

- Stores audit trail of changes made to a revalidated program
- Each log-in attempt is timestamped and kept for viewing

Certified

- Arbour Group, an independent global regulatory solutions organization serving the life science industry, conducted the third-party review of Green Button Go Automation Scheduling Software

Compliance Software Tools

- Save Validated Program
- Signatures
- View Change Log
- View Security Log
- View Signature Log

21 CFR Part 11 Compliance



Green Button Go Automation Scheduling Software


Green Button Go Automation Scheduling Software is a device-agnostic software that integrates equipment into a cohesive automation ecosystem. Whether you are integrating scientific instruments or automating manufacturing processes, Green Button Go software can automate your equipment to run 24/7, increasing output, in less time, with robust data and in-depth recordkeeping.

Equipment is easily integrated into the workflow by selecting the appropriate device driver from the driver library and positioning it in the process flow using a drag-and-drop function. It's that simple!

- Workflow design and modification
- Automation simulation
- Schedule instrument usage
- Dynamic or static scheduling
- Error recovery

Hands-on Training to Ensure Success

Hands-on training for Green Button Go software is available regularly at Biosero headquarters in San Diego, California. Classes range from Basic Training to guide software installation and set-up, to Advanced Training that teaches you to modify layouts, develop complex scenarios and manage error scenarios. Developer training is also available for people who want to learn advanced scripting and to develop their own drivers and plug-ins for Green Button Go software. **For more information, visit www.biosero.com/company/training**

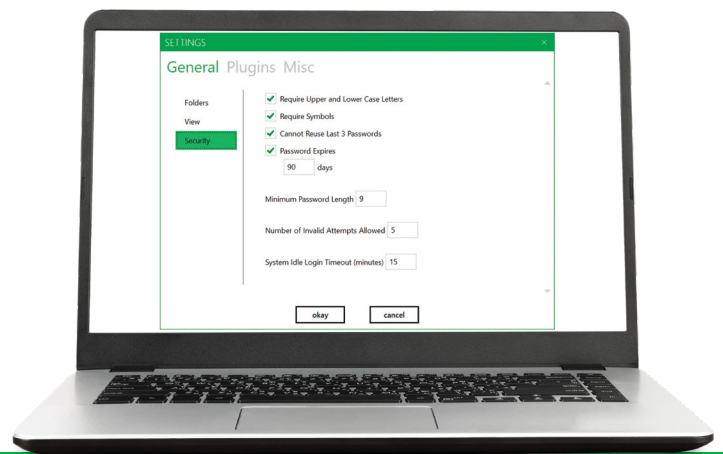

To schedule a demonstration of the Green Button Go software please call (858) 880-7376 or email info@biosero.com
www.Biosero.com/GreenButtonGo

 **BIOSERO**