

# Daniels Sharpsmart

The world's safest and most  
environmentally responsible  
**Sharps Management System**

CLINICALLY PROVEN



Safer. Simpler. Smarter.  
Greener.



# Sharpsmart Containers

## REGULAR



### S14

14 Litre Regular Sharpsmart. Ideal for Patient Room use and lighter sharps volume areas

#### Fill Capacity: 6.5 L

280mm (H) x 395mm (W) x 210mm (D)



### S22

22 Litre Regular Sharpsmart. Best suited to moderate volume areas

#### Fill Capacity: 14.5L

440mm (H) x 395mm (W) x 210mm (D)



### S32

32 Litre Regular Sharpsmart. Developed for high volume areas

#### Fill Capacity: 23.5L

610mm (H) x 395mm (W) x 210mm (D)

## ACCESS PLUS



Features a full opening lid without a tray to provide full access to the container

### S22 AP

22 Litre Access Plus Sharpsmart. Designed for secure areas where unrestricted deposit of long or bulky sharps is required

#### Fill Capacity: 14.5L

440mm (H) x 395mm (W) x 210mm (D)



### S32 AP

32 Litre Access Plus Sharpsmart. Designed for secure areas such as operating rooms with high volumes of long or bulky sharps and instruments.

#### Fill Capacity: 23.5L

610mm (H) x 395mm (W) x 210mm (D)

**According to research, a 200 bed hospital using Daniels Sharpsmart can save 1.5 tonnes of plastic waste to landfill annually.**

Compass Medical Waste Services joined forces with Daniels Health to bring a first-world, custom-made, reusable sharps solution to Africa. The Daniels Sharpsmart container solution promises safety, reduced environmental impact, cost efficiency and streamlined processes for the healthcare industry.



# Sharpsmart Accessories


## Cartsmart

Enables convenient, mobile mounting of Sharpsmart containers.




## Static Floor Stand

All Sharpsmart Regular and Access Plus Sizes.

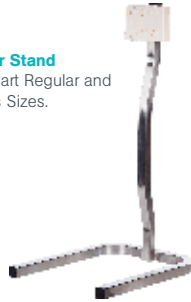

## Bench Bracket - S14 and S22

Provides secure stabilisation of S14 and S22 Sharpsmart containers on a benchtop.

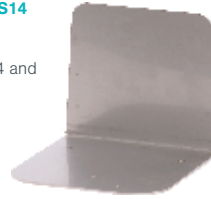

## Locking Wall Bracket

(Left Hand, Right Hand)  
For use in areas where secure, key lock mounting of Sharpsmart containers is required.



## Access Plus Wall Bracket

Designed for wall mounting of the Sharpsmart Access Plus range to allow 90-degree lid opening.



## Standard Wall Bracket

Universal bracket enabling wall mounting of Sharpsmart range. Bracket kit comes with screws and hole-plugs.




## Mounted Cradle

Identical concept and design as Standard Cradle, but with additional mounting plate to attach Cradle to a wide variety of trolley designs.



## Standard Cradle

Easy to install design simply hooks over rail on trolley to allow Sharpsmart containers to face inwards or outwards.

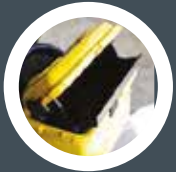

## WIDE OPENING

The optimal access geometry of the lid opening with a wide aperture. Minimises disposal-related injuries & restricts hand access to the contained waste.




## CHECK CONTENTS LEVEL

A clearview window identifies the fill level of the container, a feature unique to Sharpsmart containers.




## ELIMINATES OVERFILLING

Once the container is full, the in-built safety tray is activated into an upright closed position to restrict access and overfilling.

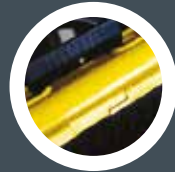

## MAXIMUM STRENGTH HARDENED PLASTIC

The medical-grade hardened plastic casing of the Sharpsmart container is completely impenetrable by contained sharps. The Sharpsmart container has been proven to withstand an amazing 500 cycles, a testament to its durability, sterility and re-usability.

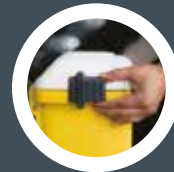

## CANNOT BE REOPENED

Side locks are tamper-proof when engaged. This safety feature prevents unauthorised access to the container contents.



## NO LEAKS

A liquid seal fixed around the rim of the container lid eliminates leakage of fluids during transport. Container leakage is a common hazard with sharps containers that aren't equipped with a leakproof seal.

## Safe disposal in a single action.

The Sharpsmart container has been strategically engineered to minimise risk by enabling single action sharps disposal:

- ▶ The maximum sensitivity safety tray is activated by the weight of the waste with a single swing action.

- ▶ The safety-tilt mechanism releases the waste safely into the container away from the user.

- ▶ When static or in motion, the swing-action safety tray effectively conceals the aperture of the container to ensure no hands ever come into contact with contaminated sharps.

- ▶ The gravity balanced mechanism of the tray resets to an open safe-disposal position after each use until the container is full.



# Across its 20+ year global history,

the Daniels Sharpsmart container has radically changed the acceptable standard of 'Safe Sharps Disposal'.

Protecting healthcare workers from needlestick injuries and diverting millions of kilograms of plastic from landfill every year, we are proud ambassadors of a product that quite literally Saves Lives.

## OUR 20 YEAR GLOBAL IMPACT



Prevented an estimated  
**85,000 needlestick injuries**



**98,732,859 single-use sharps  
containers** diverted from landfill

Safer. Simpler. Smarter.  
**Greener.**

For a free, no obligation demonstration contact us on  
**031 267 9700** or email **[daniels@compass.za.net](mailto:daniels@compass.za.net)**



[danielscompass](#)



Daniels Sharpsmart Compass