


NEURO HiTRONIC— Electrohydraulic Stance Phase Control

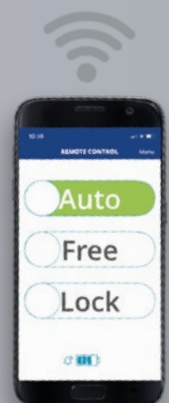

NEURO HiTRONIC – The Orthosis Control

Intelligent Motion Sensors

The controller of the orthosis has motion sensors that detect the lower leg's movements and can tell if the user is standing or walking. The control immediately reacts to changes in walking speed and step length.

Lightweight, Flat and Compact

The new controller is divided into two parts via a jointed connection and can thus be mounted anatomically and discreetly in the lower leg shell.

Various Operating Options

Using the User app on the smartphone (iOS/Android), Apple Watch or Bluetooth remote control



NEURO HiTRONIC – Electrohydraulic Stance Phase Control

Advantages for the Patient

- Infinitely locking. The system knee joint locks at any knee flexion angle and thus enables fluent motion sequences.
- Silent locking and unlocking.
- Compact and lightweight. The system knee joint can be comfortably worn under the trousers and offers a high level of wearing comfort.
- Natural gait. Slow or fast, but most importantly safe and fluent walking even when alternating between short and long steps. This gives a new walking feeling and confidence in the orthosis.
- Easy controlling. The NEURO HiTRONIC can be operated with the remote control, smart-phone (iOS or Android) or Apple Watch.
- Permanently locked. The system joint can be permanently locked in case of increased safety requirements. The NEURO HiTRONIC can also be permanently unlocked for activities (e.g. cycling).
- High performance rechargeable battery. An orthosis with the NEURO HiTRONIC system knee joint can be used all day without any problems. In Free mode: 24 hours. In Auto mode: 26,000 strides. The use in Lock mode does not require power.

Can be worn comfortably under trousers

Silent locking and unlocking

Infinite, permanent stance phase control

Intelligent FIOR & GENTZ microhydraulics

3 switch modes: Free, Lock and Auto

Rotary switch for permanent locking and unlocking


Safety thanks to Microhydraulics

The NEURO HiTRONIC is the first electrohydraulic system knee joint with stance phase control from FIOR & GENTZ. This technology, which was specially developed for orthoses, enables patients to move freely without having to think about the function of the orthosis all the time. Thanks to the intelligent motion sensors and innovative microhydraulics, the NEURO HiTRONIC supports a natural gait, thus avoiding incorrect posture and preventing damage to the musculoskeletal system.



Watch Function Video

Scan the QR code for more information

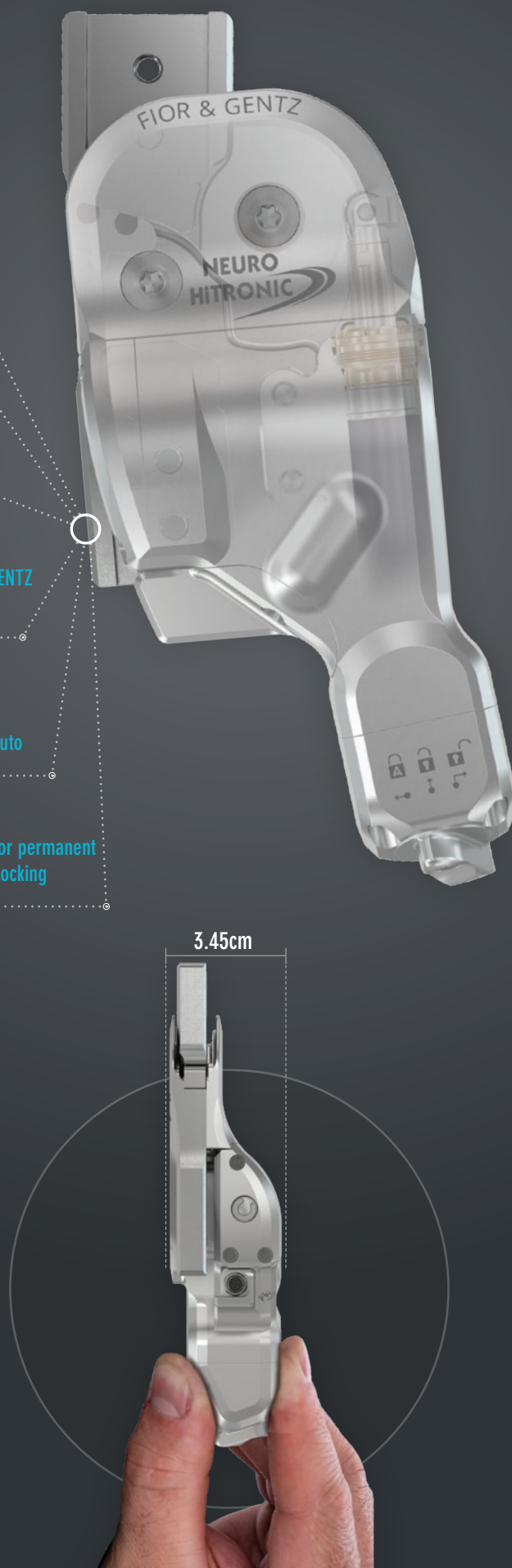

NEURO HiTRONIC
Technology
for People



Getting from A to B in everyday life is no longer a problem. The NEURO HiTRONIC offers security and flexibility in every situation.



Permanent locking and unlocking of the system joint via the rotary switch



“

Thanks to the intelligent control of the NEURO HiTRONIC, even walking on uneven terrain is easy.

”



Fits comfortably under the trousers



//

The NEURO HiTRONIC is
ideal for running around
together in the garden.
Thanks to the infinite
locking, you easily forget
it's there.

//

NEURO HiTRONIC

Easy Mounting. Easy Handling.

The electrohydraulic NEURO HiTRONIC system knee joint offers orthotists unique advantages when producing and adjusting the orthosis. Regardless of the complex technology, such as the intelligent microhydraulics and the innovative motion sensors, the mounting of the system knee joint is easy. Demounting is just as easy, for example in preparation for maintenance.

Advantages for the Orthotist

- No certification required. The system knee joint can be mounted in your own workshop without additional training.
- Expert app. Extensive adjustments on an unlimited number of NEURO HiTRONIC orthoses can be made conveniently in the workshop, for example, using a tablet.
- Unilateral treatment without supporting joint possible. This enables the production of an orthosis with the best possible function and the least material usage.
- Integrated swing phase support. The joint has an integrated, non-adjustable propulsion support that supports the end-grade knee extension.
- Easy mounting in the orthosis. As always, the orthosis is produced according to the proven FIOR & GENTZ production technique.



Easy mounting and demounting of the functional unit. The NEURO HiTRONIC can be converted to the NEURO TRONIC by replacing the system joint.





**Would you like to produce an orthosis with a
NEURO HiTRONIC system knee joint for your patient?**

Use the Orthosis Configurator to independently select the necessary system components for an orthosis with the NEURO HiTRONIC.

The Orthosis Configurator determines the appropriate system components for your patient by taking the patient data and the load capacity into account.



**Orthosis
Configurator**

www.orthosis-configurator.com

FIOR & GENTZ

Gesellschaft für Entwicklung und Vertrieb
von orthopädietechnischen Systemen mbH

Dorette-von-Stern-Straße 5
21337 Lüneburg (Germany)

☎ +49 4131 24445-0
☎ +49 4131 24445-57

✉ info@fior-gentz.de
🌐 www.fior-gentz.com

FIOR & GENTZ
ORTHOPÄDIETECHNIK MIT SYSTEM