

Remain mobil - even with an injured foot

**More mobility for more
quality of life.**



MORE MOBILITY WITH SAFETY

Have you injured your foot or ankle or just had surgery?

Everyday things become a challenge. Constantly walking with crutches requires strength and certain skills. Incorrect positioning or weight-bearing will often lead to further appointments with the doctor, subsequent surgeries or a longer period of unemployability.

The ORTHOSCOOT® is a medical aid that helps you to remain safe and comfortably mobile during this period. It enables safe mobility with little effort. The injured leg is protected by the safely guided unloading and positioning in the knee support, while the patient's mobility is ensured.

- **Guided relief and ideal positioning** of the foot
- **Safe and independent mobilization**
- In the post-operative period it **actively counteracts muscle loss.**
- Relief of the **back muscles, shoulders and wrists**
- **Hands-free** while standing
- **Angle adjustment** of the medical knee pad
- The unique **crutch-holder** enables the patient to easily attach crutches to the device
- **Foldable and can be assembled and disassembled without tools.** The weight is equal to that of a small bicycle for easy transportation
- Safe, high-quality and **TÜV-tested** according to DIN ISO 11199-2 and EN-12182
- **Sustainable and climate-neutral rental from medical supply retailers**

THE ORTHOSCOOT® IN DETAIL



Steering unit as carrying handle for easy transport



Transport compartment



Ergonomic pad with Cutout for „turning on the spot“ in narrow spaces




Toolless assembly and disassembly



Reflectors all around for more visibility




Crutch-holder

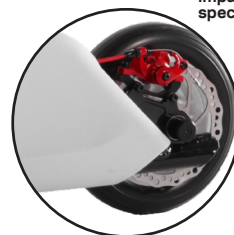

Shock- and impact-absorbing special tyres

Handbrake and locking brake

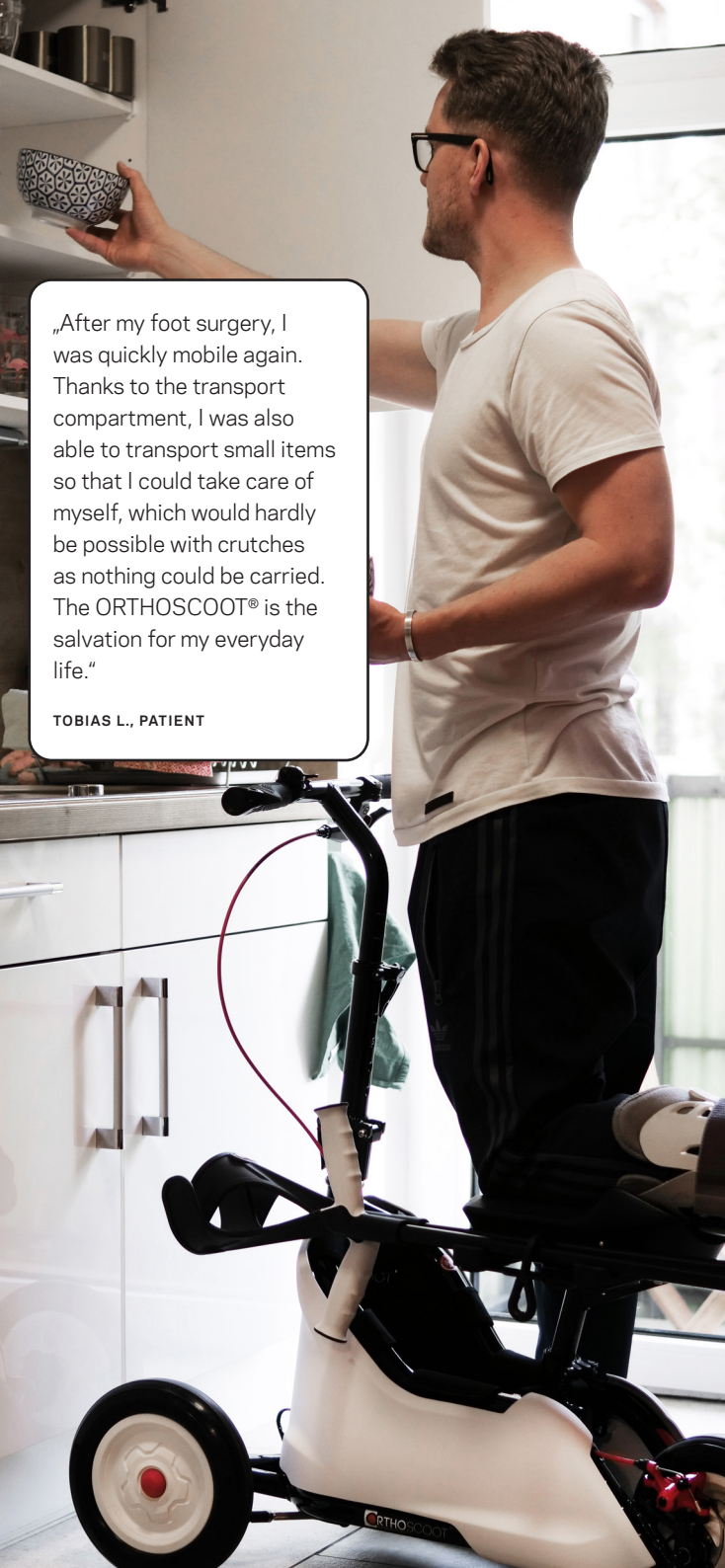A man with short brown hair and glasses, wearing a white t-shirt and black track pants, is standing on a white and black ORTHOSCOOT mobility device. He is in a kitchen, reaching up to place a small bowl with a black and white geometric pattern on a white shelf. The device has a large white front wheel and a smaller black rear wheel. A red cable is attached to the device. The background shows a kitchen counter, white cabinets, and a window with a view of greenery outside.

„After my foot surgery, I was quickly mobile again. Thanks to the transport compartment, I was also able to transport small items so that I could take care of myself, which would hardly be possible with crutches as nothing could be carried. The ORTHOSCOOT® is the salvation for my everyday life.“

TOBIAS L., PATIENT

RENT THE ORTHOSCOOT®

See your doctor

You will receive a prescription /
prescription with the appropriate
indication

1

How to rent

Go to your local medical supply
store with your prescription

2

Usage

Be mobile – enjoy everyday life
despite being injured

3

Return

Return the ORTHOSCOOT® to
the medical supply store or in the
nearest parcel store

4

Further questions? We are happy to help.

Do not hesitate to contact us! Mo-Fr from 9:00 to 17:00 (GMT+01:00) Berlin



International hotline
+49 (0) 821 710199 20



More information
orthoscoot.com

Indications - When is it suitable to use the ORTHOSCOOT®?

Limitation of mobility as a result of chronic or temporary functional limitations in the foot, ankle and/or distal lowerleg and sufficiently preserved coordination such as after/at:

- Fracture of ankle, midfoot or forefoot
- Lesion of the ankle supporting tendons
- Severe soft tissue damage of lower extremity
- Complex lesions of the ligaments around the ankle
- Amputations at the level of distal lower leg, ankle, hindfoot or toes
- Revision surgery of the lower extremity
- Injury to the lower extremity with simultaneous previous damage to the upper extremity
- Osteochondral lesions of ankle or tarsal
- After operations on the heel bone or other tarsal bones
- After significant surgery on forefoot
- In patients with overload reactions of the bone („Bone Bruise“).
- Arthrodesis of the ankle joint
- After correction of severe foot deformities like Clubfoot, flat foot, cavovarus foot

Contraindication - When is it not suitable to use ORTHOSCOOT®?

Under certain circumstances, the use of an orthopedic scooter is not suitable. Please ask your treating physician.

Technical Data

Model type	NH1 (SN starting with NH111)
Weight	14,1 kg
Length / Width / Height	920 mm / 580 mm / min. 870 mm - max. 960 mm (36,2 in / 22,8 in / min. 34,6 in - max. 37,8 in)
Range of use	Body height 140-200 cm (55,1 in-78,7 in), Body weight max. 150 kg
Noise Level	20-30 decibel on a smooth surface
Packaged dimensions	L 820 mm / W 470 mm / H 510 mm (L 32,3 in / W 18,5 in / H 20,1 in)
Girth	2780 mm (109,4 in)

Accessories

ORTHOSCOOT® lock, ORTHOSCOOT® 4th-wheel-kit

Transport and shipping

The ORTHOSCOOT® can be folded for transport in cars, airplanes and trains. For shipping, we refer exclusively to the ORTHOSCOOT® shipping box.

Maintenance and cleaning

The ORTHOSCOOT® is to be cleaned and maintained at regular intervals. Basic cleaning and disinfection are to be carried out prior to each change of user. Maintenance and cleaning may only be executed by the manufacturer. ORTHOSCOOT® may not be cleaned with aggressive, corrosive or abrasive substances. Do not use steam cleaners, hard and/or abrasive brushes or sponges. For disinfection, we refer to approved, commercially available products, but only those containing no chloride additives.

MDR (EU) 2017/745

DIN EN ISO 13485

ISO 9001

Aid number requested



Property rights

German patent application: 1020181166125

German utility model: 2020181039342, 2020171021378, 2020151064241, 2020151064233

US patents: US 10973729, US 9962312

EU patents: 3 173 059

Design protection: WO-DM/083874, US-D 770 944

ORTHOSCOOT®



ORTHOSCOOT GmbH

Gessertshausener Straße 2

86356 Neusäß

T +49 (0) 821 71 01 99-20

M info@orthoscoot.com



Climate Neutral
Company

orthoscoot.com



MIT SICHERHEIT
MEHR MOBILITÄT