



ENFit[®] PHARMACY SOLUTIONS

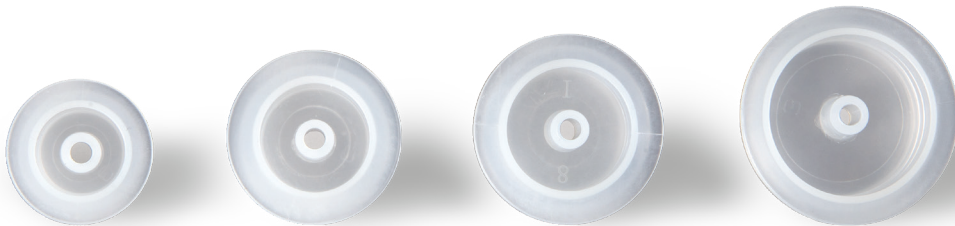

Pharmacy ENFit® Syringes

- Pharmacy syringes available in clear and amber.
- Convenient sizes of 1 mL to 60 mL in clear and amber syringes with easy-read metric scale gradients for accurate, precise dosing.
- Clinically accepted for the delivery of standard liquids.
- Complies with FDA 510K clearance in the delivery of medication.
- Comes with self-righting caps to support work flow. Capping the syringe with the self-righting caps enhances secure handling and transporting of filled syringes.
- Amber pharmacy ENFit® syringes comply with USP 40-NF 35 <67 1>, utilization Test Parameter for Spectral Transmission. USP 40-NF is the testing methodology that governs light filtration and attenuation for amber-colored containers, including syringes. USP 40-NF 35 specifies test parameters for the performance criteria of observed light transmission for plastic containers, including amber syringes.

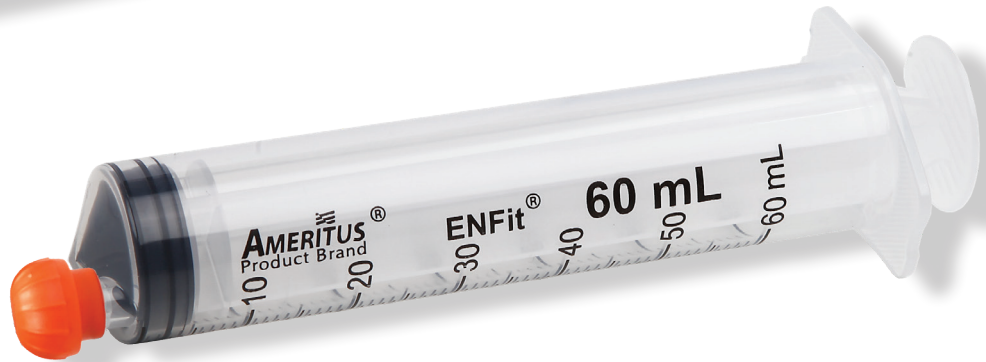

Syringe Sizes Available as Shown Above:
1mL, 3mL, 5mL, 10mL, 20mL, 30mL, and 60mL

ENFit® Tip Configurations

- Follows the GEDSA Position Statement (May2015) that recommends the use of enteral devices with ENFit® connectors to eliminate device misconnections.
- Complies with the ISO 80369-3, specific to enteral devices, with specific tip measurements to reduce the likelihood of tubing misconnections.
- Low Dose Tip. 1 mL, 3 mL, and 5 mL are designed with the Low Dose Tip. The Low Dose Tip "adds an internal male lumen to the standard ENFit® female syringe"¹. "This mimics the functionality of traditional male oral/enteral syringe designs."² Thus the accuracy of the LDT design is substantially equivalent to standard orientation (male) enteral/oral syringes.³



^{1,2,3} GEDSA, "Reducing the Risk of Medical Device Tubing Misconnections: ENFit® Low Dose Tip Syringe Review," Q2 201

Bulk Packaged Ameritus ENFit® Pharmacy Syringes			
Amber Barrel with White Gradient Markings			
Product Code	Description	UOM	QTY
ENF1-KS-LDBA	1mL*	CA	500
ENF3-KS-LDBA	3mL*	CA	500
ENF5-KS-LDBA	5mL*	CA	500
ENF10-KSBA	10mL	CA	500
ENF20-KSBA	20mL	CA	200
ENF30-KSBA	30mL	CA	200
ENF60-KSBA	60mL	CA	200

ENFit® is a federally registered trademark of GEDSA, used with permission.
Ameritus® Product Brand is a registered trademark of Kentec Medical, Inc.

Bulk Packaged Ameritus ENFit® Pharmacy Syringes			
Clear Barrel with Black Gradient Markings			
Product Code	Description	UOM	QTY
ENF1-KS-LDBC	1mL*	CA	500
ENF3-KS-LDBC	3mL*	CA	500
ENF5-KS-LDBC	5mL*	CA	500
ENF10-KSBC	10mL	CA	500
ENF20-KSBC	20mL	CA	200
ENF30-KSBC	30mL	CA	200
ENF60-KSBC	60mL	CA	200

*Denotes Low Dose Tip Syringe Design

Press In Bottle Adapters for ENFit® Connectors



Press In Bottle Adapters for ENFit® Pharmacy Syringes			
Individual Blister Packs			
Product Code	Description	UOM	QTY
AMP-PIBA-1	12-15mm Outside Diameter	BX	100
AMP-PIBA-2	14-16 Outside Diameter	BX	100
AMP-PIBA-3	17-20 Outside Diameter	BX	100
AMP-PIBA-3.5	19-23 Outside Diameter	BX	100
AMP-PIBA-4	24-26 Outside Diameter	BX	100

The Ameritus Press In Bottle Adapters make any bottle ready to dispense. Insert into bottle, attach your ENFit syringe, and draw your dose of medication. We stock four sizes that fit most standard bottles. Once the medication is drawn the bottle caps can be placed back on the bottle with the adaptor. Single packed. Single use.

Distributed by

Kentec Medical, Inc.

17871 Fitch Irvine, CA (949) 863-0810 FAX (949) 833-9730

TOLL FREE - USA Only (800) 825-5996

sales@kentecmedical.com | www.kentecmedical.com



Visit us on Facebook, YouTube and LinkedIn!