

The UK's Leading **Abduction Prevention** Solution.

Xtag Medical are the UK/NHS' leading supplier of baby tagging security systems. They are prioritised by over 60 NHS trusts for their affordability, ease of use and long-term reliability in fast-paced neonatal environments.

Our technology provides the highest level of security for new-borns and our systems are designed and installed by patient safety specialists based in the UK to meet each hospital's unique requirements.

————> [Watch Solution Video](#)



Solution Overview.....1

Contents.....2

Solution Benefits.....3

Baby Tags.....4

XTM6™ Baby Tags.....5

Tag Features.....6

Resources.....13

Patient Protection.....7

Security Triggers.....8

Abduction Prevention.....9

Patient Management.....10

Management Dashboard.....11

Dashboard Functions.....12

Solution Benefits



Award-winning Security

Award-winning child safety technology ensures the greatest protection for new-borns.



Approved by Regulators

Recommended by external regulators such as the Care Quality Commission (CQC).



Scalable

Fully scalable across multiple floors and hospital buildings.



Designed for Midwives

Designed to provide maximum practicality and ease of use for maternity staff.

Safeguarding Over **300,000 Babies** in over **60 NHS Trusts** Annually.

Essential in modern day healthcare, XTM6™ Electronic Baby Tags provide the industry's **greatest protection from the risk of abduction**, ensuring peace of mind for families and midwives.



XTM6™ Baby Tags

Real-time Protection

Baby tags transmit their location and safety status **every second**, providing **real-time intelligence** on a new-born's safety throughout the day.

IP Ratings	IP67, IP65
Enclosure	White
Dimensions	39 x 30 x 10mm (L x W x H)
Weight	12 Grams

Contact info@xtag.co.uk
for a Free Sample.



Tag Benefits



Comfortable

Small, light and round for a snug fit on a baby's ankle.



Safe

No sharp edges and extremely durable material.



Long Life

Battery lasts up to three years without the need to charge.



Easy Clean

Easy to clean and sterilise ready for re-use.



Water-resistant

IP67 water-resistant, so no need to remove for bathing.



Hypoallergenic

Our tags and straps won't cause allergic reactions.


Award-Winning **Patient Protection** Technology.

Award-winning child safety technology detects when a new-born's safety is at risk. Our systems are designed to integrate with **existing security infrastructure within the hospital**, including Access Control, CCTV and elevator systems.



Security Triggers

A ward-wide sounder alarm occurs when there is a serious threat to a new-born's safety.

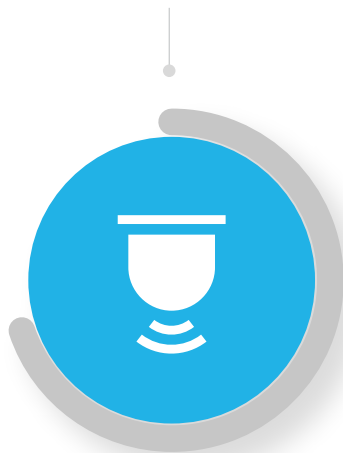


1

Abduction Alarm

2

Monitor Alert

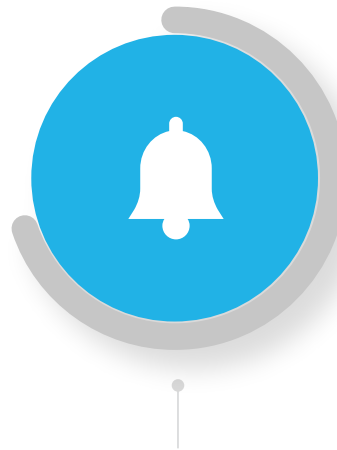

An alert is sent to the touchscreen monitor for midwives/security to investigate a safety issue.















The system forces local exits to lock, preventing a potential abductor from leaving the ward.



3

Exit Lockdown

Abduction Protection

Protection Function	Description	Security Triggers
Exit Protection	Occurs when a baby approaches an <u>open</u> exit.	  
Loiter Protection	Occurs when a baby approaches a <u>closed</u> exit.	 
Tag Off Protection	Occurs when a tag is removed without cutting the strap.	  
Lost Signal Protection	Occurs when a tag is not seen by the system for a few seconds.	  
Tamper Protection	Occurs when a strap is cut before a patient has been discharged.	  

All-In-One **Patient Management** Dashboard.


The Touchscreen Client is an **all-in-one patient management system** that enables complete patient management on a single, easy to use dashboard. Comes with **barcode scanner** for quick enrolment of patient tags.



Management Dashboard

Alarm status alerts when a baby's safety is at risk

Last seen location status updated every second

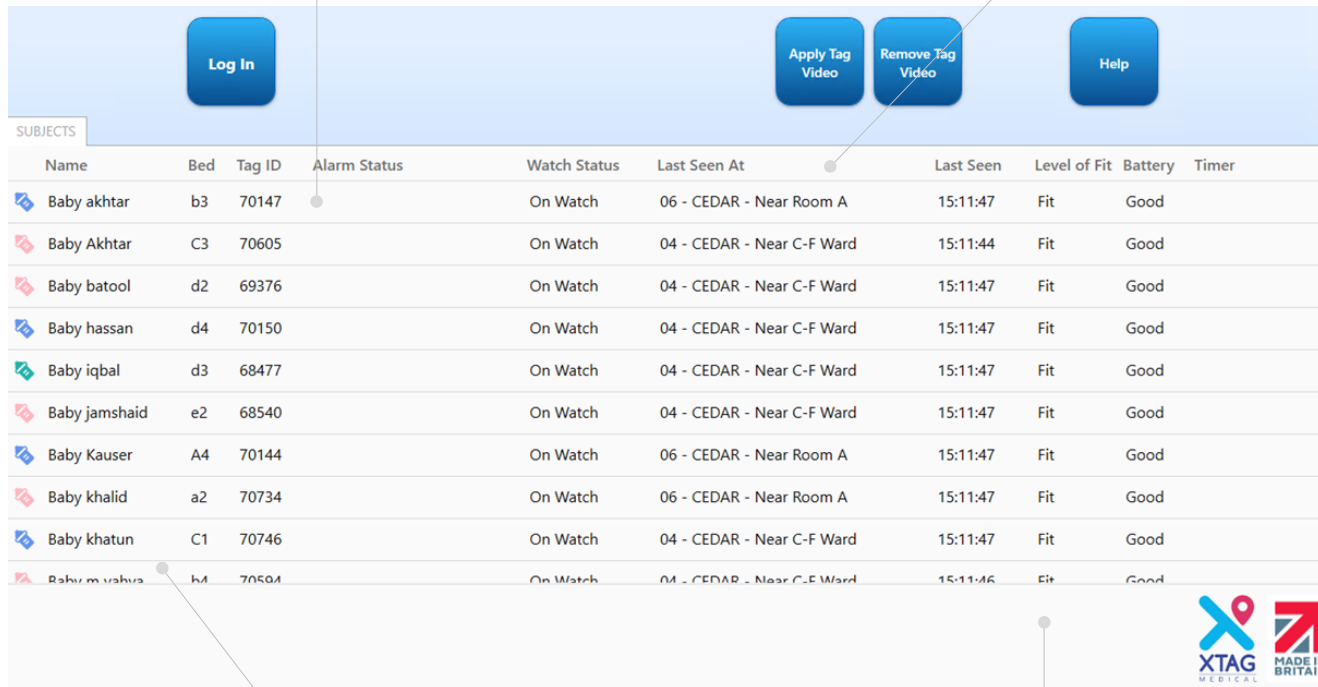

The dashboard interface includes a top navigation bar with buttons for 'Log In', 'Apply Tag Video', 'Remove Tag Video', and 'Help'. Below this is a table titled 'SUBJECTS' with the following columns: Name, Bed, Tag ID, Alarm Status, Watch Status, Last Seen At, Last Seen, Level of Fit, Battery, and Timer. The table lists ten babies, all with a status of 'On Watch' and 'Fit' level of fit. The 'Last Seen' column shows timestamps like '15:11:47'.

Each baby given a customisable unique ID

Ensure tags are worn correctly with level of fit indicator

Dashboard Functions



Alarm Status

Shows patient location and safety status in real time.



Transfer Function

Enables tagged baby transfers between wards.



Off Watch Function

Enables transfers to examination rooms, treatment rooms etc.



Discharge Function

Remove discharged patients from the system.



Valuable Data

View audit logs, analytics and archive valuable data.



Secure Access

Unique access for each user and advanced supervisor access.

Resources

For more information on Xtag Medical visit: <https://xtag.co.uk/>

For more information on our solutions visit: <https://xtag.co.uk/Solutions/>

Request a live demonstration or a free quote: <https://www.xtag.co.uk/contact/>

Get in touch:

Email: info@xtag.co.uk

Phone: [0333 340 2390](tel:03333402390)