


Above Tarsus (Chevy) Prosthetic Measurement Form

Circle: Left/Right

Pet Name: _____ Breed: _____ Age _____ Lbs: _____

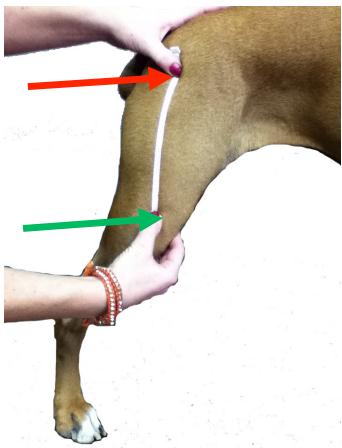

1. Hip to Stifle _____

greater trochanter of femur
apex of fibular head




4. Thigh Circumference

R _____ L _____



2. Groin to Stifle _____


* femoral pulse*
medial condyle of femur

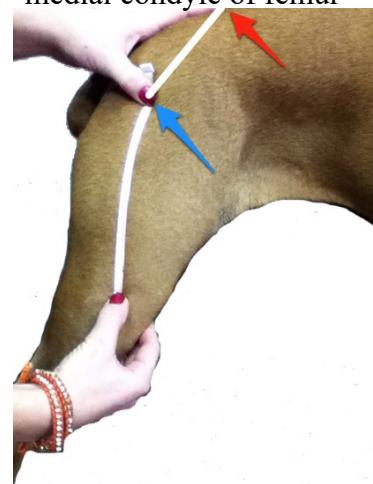

5. Greater trochanter to LS
junction

R _____ L _____

* Greater Trochanter*
Lumbosacral Junction

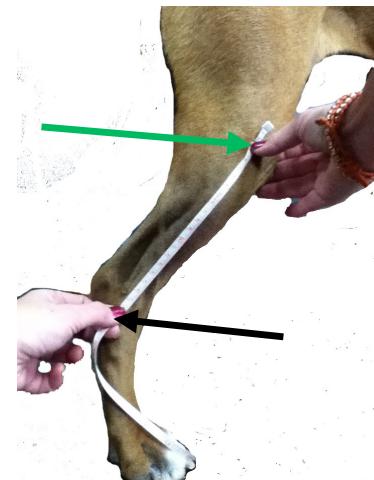

3. Stifle to end of limb _____


apex of fibular head
end of residual limb



6. Greater Tibial Tuberosity
to Top of Patella

7. End of limb to
ground _____

ORTHOPETS®

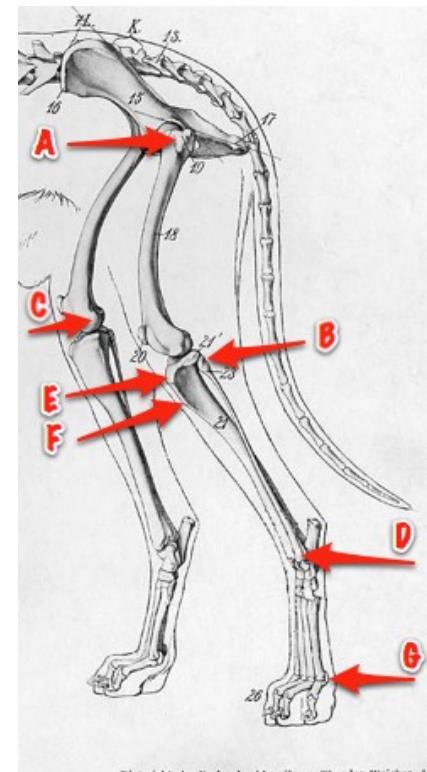

A: Greater Trochanter of Femur

B: Fibular Head

C: Medial Condyle of Femur

D: Lateral Malleolus

E; Greater Tibial Tuberosity

F: Lesser Tibial Tuberosity

G: 5th Metatarsal Head

7. Standing angle of
tarsus of "Good" leg
