

Neo Cot

Bassinet



Neo Cot

A bassinet to bring mother and baby closer together.

A bassinet to encourage interaction between mother and baby. Developing mother/baby bonding by creating a transparent cot which promotes eye contact and visibility of the baby.



* Shown here is a High-Low stand.

“Safety”

The transparent cot can be locked, minimizing the risk of the baby falling.

The transparent cot can be locked at 4 points underneath the mattress to ensure baby safety.



“Interaction”

The height of the HL stand can be adjusted to promote eye contact between mother and baby.

The lowest position of the HL stand is 62cm which enables the cot to be adjusted to bedside height, allowing bonding to take place between mother and baby.

The Neo Cot has been designed to support the caring model of both mothers and nursing staff, paying special attention to safety and interaction.

Caster-lock from any direction

Four-wheel locking mechanism ensures the cot remains stable when interacting with the baby.



Mattress can be tilted using a simple lever situated at the front of the cot

Stepless smooth tilting is possible while watching a baby.



Lightweight and easy to move

Large dual wheel casters make it easy to operate and move the cot.



A combination of options can be selected to meet a variety of nursery requirements.

Enhances operation efficiency for nursing staff

Variety of peripheral devices can be attached to the IV Pole (option), which enables a staff to operate the devices immediately in a large work space.

Drawer

Pulse oximeter can be put on the slide tray which is suitable for practical use.



Sliding open the top slide tray



Opening the drawer

IV Pole

IV Pole can be attached to each corner.

Basket

Small items such as disinfectant and thermometers can be stored.



Carefully selected material ensures improved ease of cleaning.

Increased ease of cleaning of the product

The main body of the Neo Cot is made of resin which allows for a variety of different cleaning materials to be used.



The mattress is washable

Highly breathable and washable material.



Specifications

		Neo Cot HL Stand	Neo Cot Fixed Stand
Atom's Code		26408	26409
Tilting Angle		0° - 12°	
Mattress	Material	Polyester-elastomer (BREATHAIR®)	
	Method of disinfection	Dipping (avoid using dry heat sterilization)	
Caster		Dia. 10 cm (Dual caster w/stopper)	
Load Capacity		10 kg (on mattress)	
Weight		Approx. 27 kg	Approx. 25 kg
Dimensions		W960 x D450 x H620 - 820 mm	W960 x D450 x H750 mm

Options & Supplies



Atom's code: 26113

I.V Pole for Neo Cot



Atom's code: 26114

Drawer for Neo Cot



Atom's code: 26117

Basket for Neo Cot



Atom's code: 26411

Cylinder Rack for Neo Cot



Atom's code: 26123

Shade for Neo Cot



Atom's code: 26118

Mattress for Neo Cot



Atom's code: 26122

Mattress Sheet for Neo Cot (5 pcs/pkg)


Atom's code: 26112 **Transparent Cot for Neo Cot**

Atom's code: 60337 **Neo Sheets (100 pcs/box)**

Atom's code: 67100 **Baby card , blue (50 pcs/pkg)**

Atom's code: 67101 **Baby card , pink (50 pcs/pkg)**

Dimensional drawing
Unit:mm

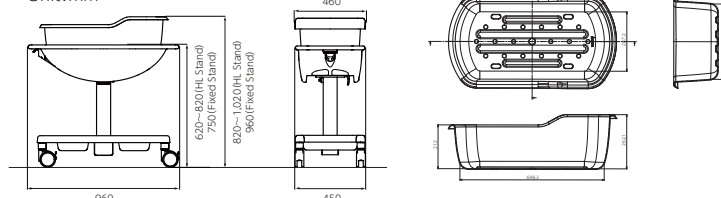

Atom Medical products are produced under strict production control in Japan.
We provide value and trust of the highest quality for the tiny baby lives in the world.

Manufacturer: ATOM MEDICAL CORPORATION

18-15, Hongo 3-chome, Bunkyo-ku, Tokyo, 113-0033 JAPAN
<http://www.atomed.co.jp>

Exporter: ATOM MEDICAL INTERNATIONAL, INC.

Iwakata Bldg., 18-16, Hongo 3-chome, Bunkyo-ku, Tokyo, 113-0033 JAPAN
Phone: + 81 3 3815 2941 Fax: + 81 3 3812 9670/4080